


Panamatic Mini Bedpan Washer Disinfector

Redefining Infection Control





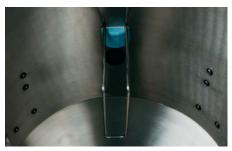







Redefining Infection Control

Our range of DDC Dolphin bedpan washer disinfectors have been developed over decades to provide you with a fast and safe human waste disposal system with passive safety measures to minimise the risk of infection. Designed to empty, wash and disinfect human waste containers such as bedpans, commode pots and urine bottles, our bedpan washer disinfectors incorporate the latest technological advances.

The Panamatic Mini bedpan washer disinfector balances size and function, whilst still offering our cutting-edge cleaning technology.

Designed for low-use environments and capable of disinfecting up to 2 reusable utensils per cycle, it allows for high-pressure cleaning in a cycle time of just 6 minutes to help maximise efficiency.

Benefits

- Reduce cross contamination of infectious materials and minimise the risk of ward outbreaks, keeping staff and patients safe.
- Incorporates the latest technological advances and able to achieve thermal disinfection levels in line with HSE periodic testing guidelines.
- Smaller capacity model to suit lower usage healthcare environments.
- Extremely thorough wash performance, eliminating 99.99% of known bacteria.



+353 (1) 443 3390



Panamatic Mini Features

Automatic scale management system

- suitable for soft and hard water areas, to prevent blockages.

Smaller capacity to suit low usage environments can clean 1 bedpan and 1 urine bottle per cycle.





Completely hands free, no touch operation

- reducing the risk of infections spreading to clinicians and further patients.



for a simple loading and emptying process.



High quality materials and components

- including antimicrobial coatings and smooth, crevice-free surfaces, reducing the presence of trapped harmful bacteria.



Antimicrobial surfaces

- to minimise cross infection and reduce the possibility of microbial growth.



Unique multi-utensil racking design

- extremely convenient with one rack to fit all standard bedpan sizes.



Rotating wash nozzles and strategically sited spray jets

- ensure a superior wash performance.



Optional programmable audio

- in your own language to ensure correct and proper use by staff.



Compact – to save space, making it suitable for all room types and sizes.



PANAMATIC

Technical Specifications

General

Stainless Steel grade: 30

 $\begin{array}{ll} \mbox{Dimensions w x h x d (lid closed):} & 500 \times 1042 \pm 3 \times 600 \mbox{mm} \\ \mbox{Dimensions w x h x d (lid open):} & 500 \times 1500 \pm 3 \times 600 \mbox{mm} \\ \end{array}$

Weight (Machine Only):70kgWeight (Packed):85kgWeight (Export):90kg

Capacity: 1 pot + 1 urine bottle

Operation cycle time: 6 minutes
Decibel rating: 49.4dBA

Electrical

110V Voltage: N/A

230V Voltage: 230V / 50Hz with dedicated 13A

breaker (Minimum)

230 Voltage: 230V / 60Hz with dedicated 13A

breaker (Minimum)

Circuit protection (230 Voltage): 30mA RCD, GFI, IEC60364 compliant Cable connection (230 Voltage): 3 core sheathed x 1.5mm². 2m long

Cable connection (230 Voltage): 3 core sheathed x 1.5mm². 2m kg 3 Phase Voltage: N/A

Circuit protection (3 Phase Voltage):

Cable connection (3 Phase Voltage):

N/A

Rated power input:

Electrical consumption per cycle:

Steam generator power:

Motor power:

Pump power:

N/A

0.425kW

Water

Water connection: 3/4* BSP male fittings required

Water pressure: (washing machine hose)
Water pressure: Cold water 0.5-10 bar
Water flow rate (Min): 10L per minute at 1 bar

Water tank capacity: 20L Water consumption per cycle: 20.4L

Waste connection: 75mm ID (3") Pushfit


Overflow connection: 22mm

Disinfection

Recommended disinfection detergents: EcoWash
Disinfection tank volume (external): 5L
Disinfection dosing per cycle: 30mL

Dimensions

Clearance for lid

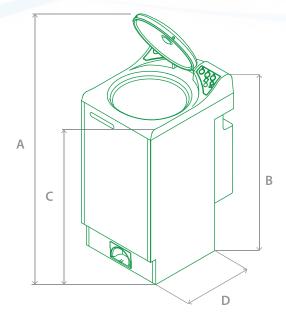

nnect	ed to	left	or r		
of t	he m	nach	ine	_	
K			I B	D	

Overflow output can be

	Millimetres			
Α	1500 ±3			
В	1042 ±3			
C	960			
D	600			
Е	766			
F	500			



^{*}Specifications subject to change







CaretuaUnit 3 Westpoint Business Park
Mulhuddart
Dublin 15
D15 R295



www.caretua.com



+353 (1) 443 3390



info@caretua.com





